Product Specification of:

Medical Oxygen Sensor

Type: M-80

Part Number: 410031

RoHS compliant by exemption, SVHC: lead
Meets the applicable requirements of ISO 80601-2-55
Comes with CE marking, further regulatory registrations available upon request
Produced under EN ISO 13485 Quality Management System

DOCUMENT PURPOSE

The purpose of this document is to present the performance specification and key features of the sensor.

This document should be used in conjunction with the Operating Manual of the instrument and the Product Safety Data Sheet of the sensor.



KEY FEATURES

Dual-cathode galvanic gas sensor that shows high signal stability at low cross sensitivities to anesthesia gases combined with superior linearity over the entire measurement range.

High Signal Stability	Superior Linearity	Wide Measurement Range	Low Cross-Sensitivities

IT Dr. Gambert GmbH .: Hinter dem Chor 21 .: 23966 Wismar .: Germany
Phone: +49 (0)3841 220 050 .: Fax: +49 (0)3841 220 052 2 .: E-Mail: info@itg-wismar.de **www.itg-wismar.de**

TECHNICAL SPECIFICATIONS

MEASUREMENT

Operating Principle:	Partial pressure electrochemical cell	
Measurement Range:	0 ... 100 Vol.% O_2	
Initial Output Signal:	11 ... 14 mV	at dry ambient air
Output Signal Difference:	1.25 ... 1.6 mV	
Response Time t90:	< 12 s	
Signal Drift (long-term):	≤ ± 1 Vol.% O_2 / month	at ambient air
Linearity Error:	≤ ± 2 %	at 100 Vol.% O_2 applied for 5 min
Zero Signal Offset:	≤ 200 µV	at 100 Vol.% N_2 applied for 5 min
Repeatability Error:	≤ ± 1 Vol.% O_2	at 100 Vol.% O_2 applied for 5 min
Influence of Gas Humidity:	- 0.03 % of signal per %rH	
Signal Temperature Compensation:	NTC on sensor PCB	
Cross-Sensitivity:	fulfills DIN EN ISO 80601-2-55	

ELECTRICAL

Electrical Connector:	3-Pin Molex
Recommended Load Resistor:	≥ 10 kOhm

MECHANICAL

Weight:	≤ 26 g
Material in Contact with Media:	PA, PPS, PTFE, Stainless Steel, NBR
Gas Connector:	fits for M16x1 DIN 13 or 5/8-24 UNEF

ENVIRONMENTAL

Operating Temperature Range:	10 ... 45 °C	
Ambient Pressure Range:	700 ... 1250 hPa	hPa
Ambient Humidity Range:	up to 100 %rH	non-condensing

LIFETIME

Expected Operating Life:	3 years	at ambient air, depending on application
Nominal Sensor Life:	750 000 Vol.%h O_2	at ambient air, depending on application

Important Note: All characteristics are based on conditions at 25 °C, 50 %rH, 1013 hPa and a gas flow of > 2.5 L_S/min. For sensor performance data under other conditions, contact ITG.

IT Dr. Gambert GmbH .: Hinter dem Chor 21 .: 23966 Wismar .: Germany
Phone: +49 (0)3841 220 050 .: Fax: +49 (0)3841 220 052 2 .: E-Mail: info@itg-wismar.de **www.itg-wismar.de**

STORAGE CONDITIONS IN UNOPENED ORIGINAL PACKAGE

Packaging:		
Ambient Temperature Range:	15 ... 25 °C -20 ... 50 °C	recommended maximum (≤ 10 h)
Ambient Pressure Range:	700 ... 1250 hPa	
Ambient Humidity Range:	50 ... 100 %rH 0 ... 30 %rH	recommended, non-condensing maximum one week
Shelf Life:		

RELATED PRODUCTS

Product	Part Number	Other Specifics
O_2 Sensor M-80	41 00 31	--

Installation in Device

A gas tight sealing is ensured when the sensor is screwed in hand-tight. Do not use any mechanical tools to install the sensor. Using excessive force may damage the sensor.

Optimal mounting position of the sensor is when the gas sensing area faces downwards. A horizontal position is acceptable. It is not recommended to use the sensor with the gas sensing area facing upward.

Connection should be made via recommended electrical and mechanical connectors only. The specified load resistance must be taken into account. Soldering the sensor will damage it and void the warranty. Please contact ITG for further information.

Since temperature has an influence on the output signal it is not advisable to place any heat or vibration sources (i.e. electric-pumps or valves, coolers, etc.) in close proximity to the cell.

Avoid proximity of the cell to any EMC radiating units as those emit frequencies which might interfere with the sensor`s electrical connection and connecting cable.

Cleaning and Disinfection

The sensor outer housing can be cleaned with a dry wipe. Do not wipe the sensor`s gas entrance side. Do not use any chemical disinfectant or sanitizer on the sensor.

Poisoning

ITG sensors are designed to operate in a wide range of environments. For optimal sensor lifetime and performance it is important that exposure to high concentrations of solvent vapors is avoided during storage, installation into instruments and operation of the sensor.

Do not use adhesives directly on or near the sensor as the solvents may cause stress corrosion on the plastic parts.

Intended Use

The electrochemical oxygen sensors for use in medical technology are used as accessories with a limited shelf life exclusively in conjunction with other medical products (ventilators, incubators, anesthesia machines, portable oxygen monitors, oxygen therapy devices and gas mixing units) so that they can fulfill the manufacturer`s intended purpose.

The oxygen sensors measure the oxygen partial pressure in gas mixtures.

Commissioning and use is carried out exclusively by expert clinic or service personnel. Direct contact with these items occurs during installation or maintenance by trained personnel.

IT Dr. Gambert GmbH .: Hinter dem Chor 21 .: 23966 Wismar .: Germany
Phone: +49 (0)3841 220 050 .: Fax: +49 (0)3841 220 052 2 .: E-Mail: info@itg-wismar.de **www.itg-wismar.de**

Stabilization Time

When installing a new sensor refer to the instrument manual for stabilization time before calibration.

If not specified otherwise wait at least 15 minutes to ensure that the sensor has stabilized in the instrument.

Calibration Interval

ITG sensors are designed to have minimal signal drift over their functional lifetime. For optimal performance and maximum measurement accuracy however they should be calibrated before each use.

If the Sensor is dropped

If a sensor is dropped, please check for visible mechanical damage of the sensor or if the sensor is leaking electrolyte. If this is the case, do not try to install the sensor into the device. Take safety precautions and immediately dispose the cell.

If the sensor shows no visible damage it should be placed in quarantine for at least 24 hours. Afterwards a follow-up check made by a two point calibration should be done.

RFI/EMI Susceptibility


ITG sensors contain metal parts and might be susceptible to RFI or EMI. Before use in MRI environments please contact ITG for further information.

Disposal

At the end of the sensors lifetime the sensor should not be disposed of in normal public waste as it may contain hazardous materials and caustic electrolyte (for more information refer to the PSDS). Please contact your local authorities for environmental legislation to relevant local waste disposal.

IT Dr. Gambert GmbH .: Hinter dem Chor 21 .: 23966 Wismar .: Germany
Phone: +49 (0)3841 220 050 .: Fax: +49 (0)3841 220 052 2 .: E-Mail: info@itg-wismar.de **www.itg-wismar.de**

Description of Pictograms

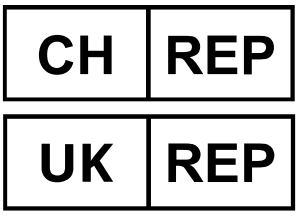
	Manufacturer
	Medical device
	Model number Customer sensor type / Customer article number
	Sensor part number of the manufacturer ITG sensor type / ITG article number
	Serial number
	Contains hazardous substances
	Indicates the temperature limits to which the medical device can be safely exposed
	Expiration Date (DD-MM-YYYY)
	Responsible representative in the territory
	Federal law restricts this device to sale by or on the order of a licensed practitioner (U.S. only)
	CE Marking with certifying notified body number

IT Dr. Gambert GmbH .: Hinter dem Chor 21 .: 23966 Wismar .: Germany
Phone: +49 (0)3841 220 050 .: Fax: +49 (0)3841 220 052 2 .: E-Mail: info@itg-wismar.de **www.itg-wismar.de**


TECHNICAL DRAWING

Sensor	ITG Part number	Color			Electrical pin assignment		
		Cap	Housing	Threaded ring	Pin 1	Pin 2	Pin 3
M-80	410031	white	white	white	plus (+)	minus (-)	plus (+)


Top view

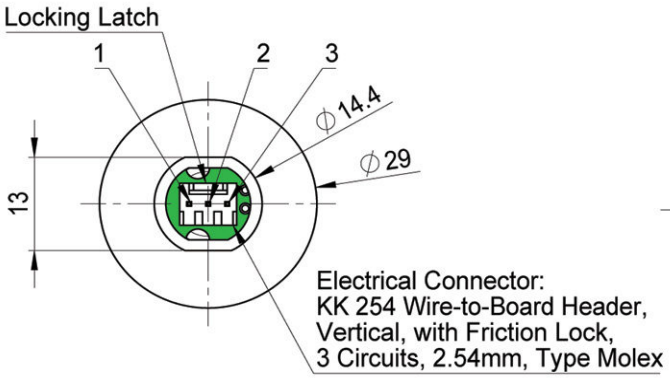

Bottom view

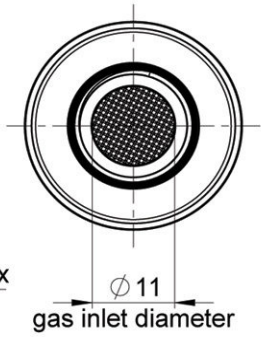

Front view

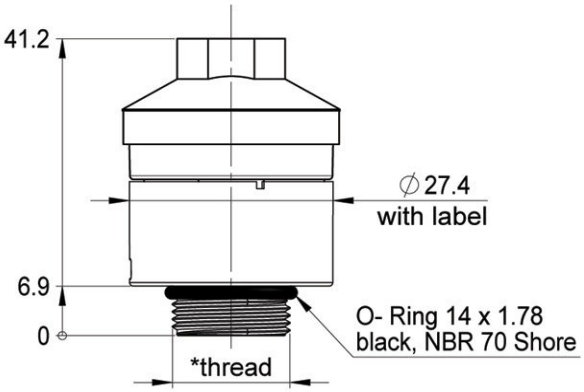

[* thread is compatible with DIN 13 M16x1 and UNEF 5/8-24]

Dimension unit: mm
Dimension tolerances: linear ±0.5 mm , diameter ±0.3 mm

IT Dr. Gambert GmbH .: Hinter dem Chor 21 .: 23966 Wismar .: Germany
Phone: +49 (0)3841 220 050 .: Fax: +49 (0)3841 220 052 2 .: E-Mail: info@itg-wismar.de **www.itg-wismar.de**